

# SYRINGE SHIELDING

New generation of personal protection for syringes

The low energy syringe shield is designed to give the best mix of personal protection and a high functionality, together with fitting the syringe perfectly. Hoy syringe protection comes in different models. The syringe shields are available with three different types of locks: Spring Lock, Twist Lock and Screw Lock.

## TUNGSTEN SYRINGE SHIELD

LEAD GLASS HIGH DENSITY (5,2 g / cm<sup>3</sup>)  
TUNGSTEN SHIELDING 2,5 mm

The tungsten syringe shield is designed for all types of syringes and for the following sizes:  
1 ml, 2 ml, 2,5 ml, 3 ml, 5 ml, 10 ml, 30 ml and 50 ml.

The scale of the syringe is fully visible through the leadglass.



## LEAD GLASS SYRINGE SHIELD

LEAD GLASS HIGH DENSITY (5,2 g / cm<sup>3</sup>)  
THICKNESS 5,5 mm.


The lead glass syringe shield is designed for all types of syringes and for the following sizes:  
1 ml, 2 ml, 2,5 ml, 3 ml, 5 ml, 10 ml, 30 ml and 50 ml.



## LOCK TYPES

### SPRING LOCK

The spring lock encloses the syringe by pushing a button to release a small stick that attaches the syringe into the sleeve.

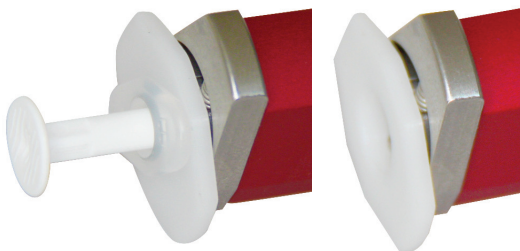

### TWIST LOCK

The twist lock encloses the syringe by twisting the syringe into the two flanges that attaches the syringe into the sleeve.

